



Technical:

1. Adaptor working voltage:AC100-240V, 50-60HZ
2. 9 settings for adjustable rotational speed, ranging from 125 to 625rpm
3. Adjustable torque, range:0.3-5.2N.cm(3-52mN.m)

Easy operating with With LCD Screen




Originals Switzerland electro-motor

More stable and energy efficient,Long service life




Top priority: patient safety, 4 different safety functions

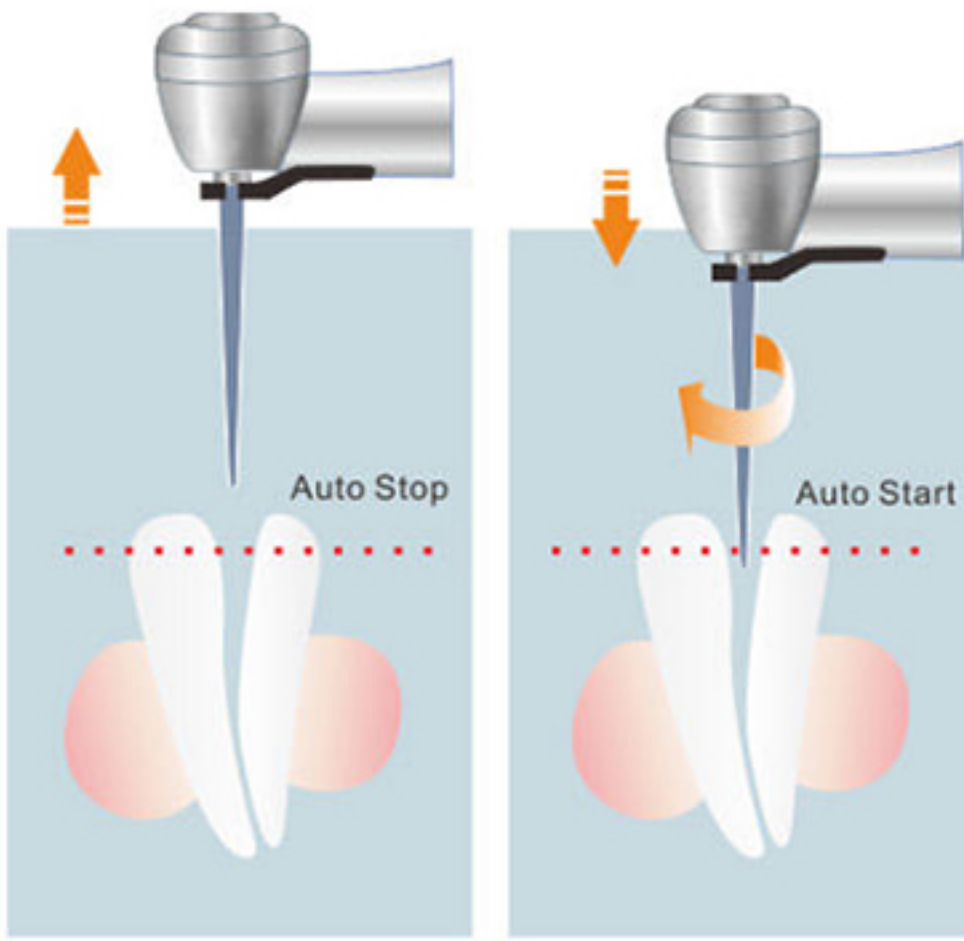


Auto start/stop

The file starts rotating auto-matically when inserted into the canal, and stops automatically upon withdrawal.

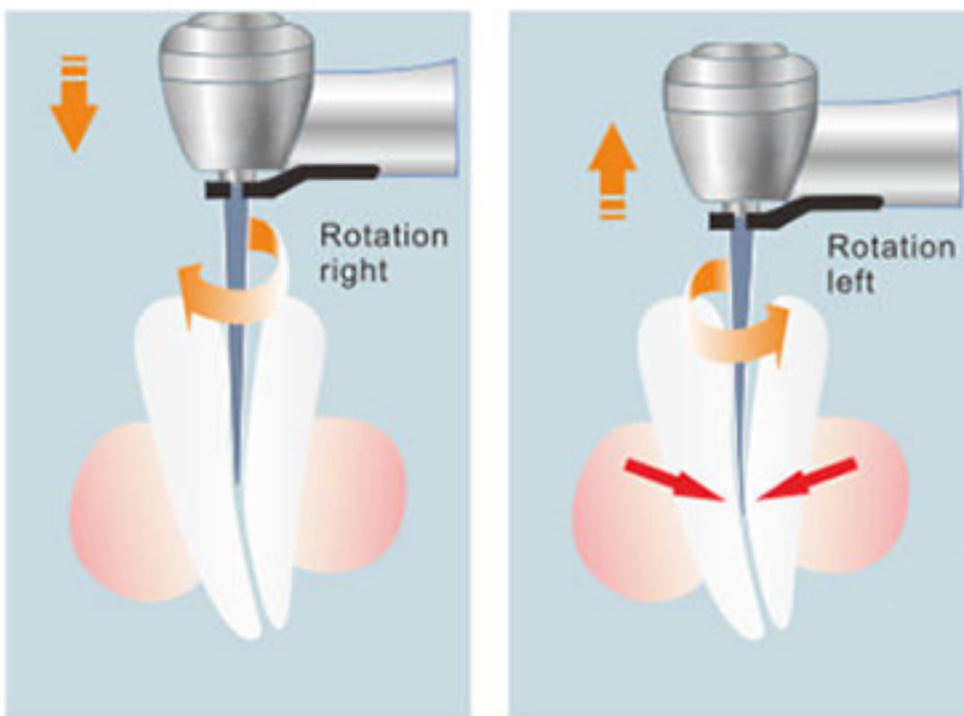


Auto torque reverse

The file stops rotating auto-matically if the torque exceeds the specified limit of one of 11 torque levels.

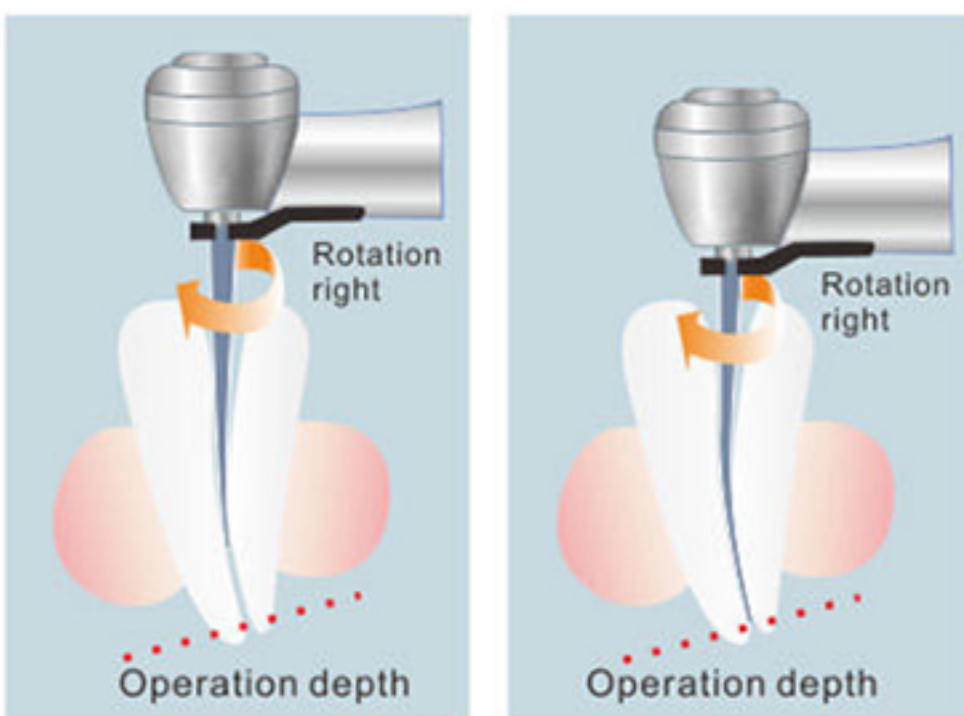

Auto apical reverse

When the file reaches at a specified measurement. The file rotation stops and reverses automatically.

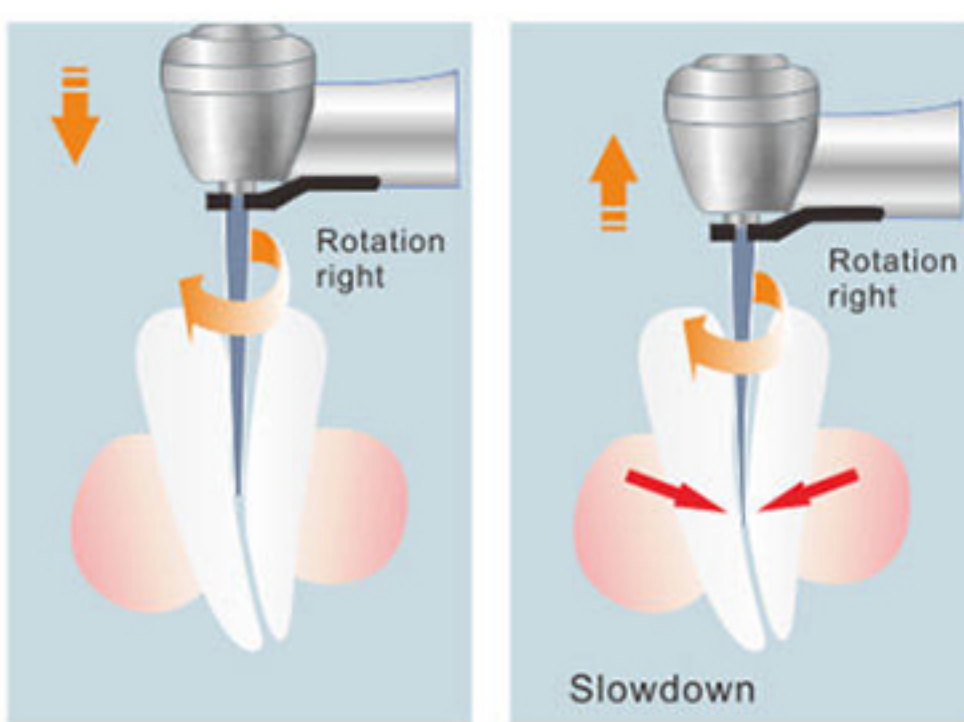

Auto apical slow down

When the file is approaching to apical foramen, the motor auto-matically slows down for safer treatment.

Feature:

1. Original Switzerland electro-motor
2. Easy operating with wide LCD screen
3. Drive-1:1 contra angle is available
4. Adjustable speed and torque
5. Auto reverse to protect the file
6. High capacity Li-ion chargeable battery and can work with battery or adaptor
7. Auto power off and memory function